


## Vascularperspectives

innovations for the interventionalist

ASAHI PERIPHERAL  
GUIDE WIRES

## ASAHI's Wire Technology Accessing the Peripheral Field

Making our superior wire technology available for superior peripheral interventions.

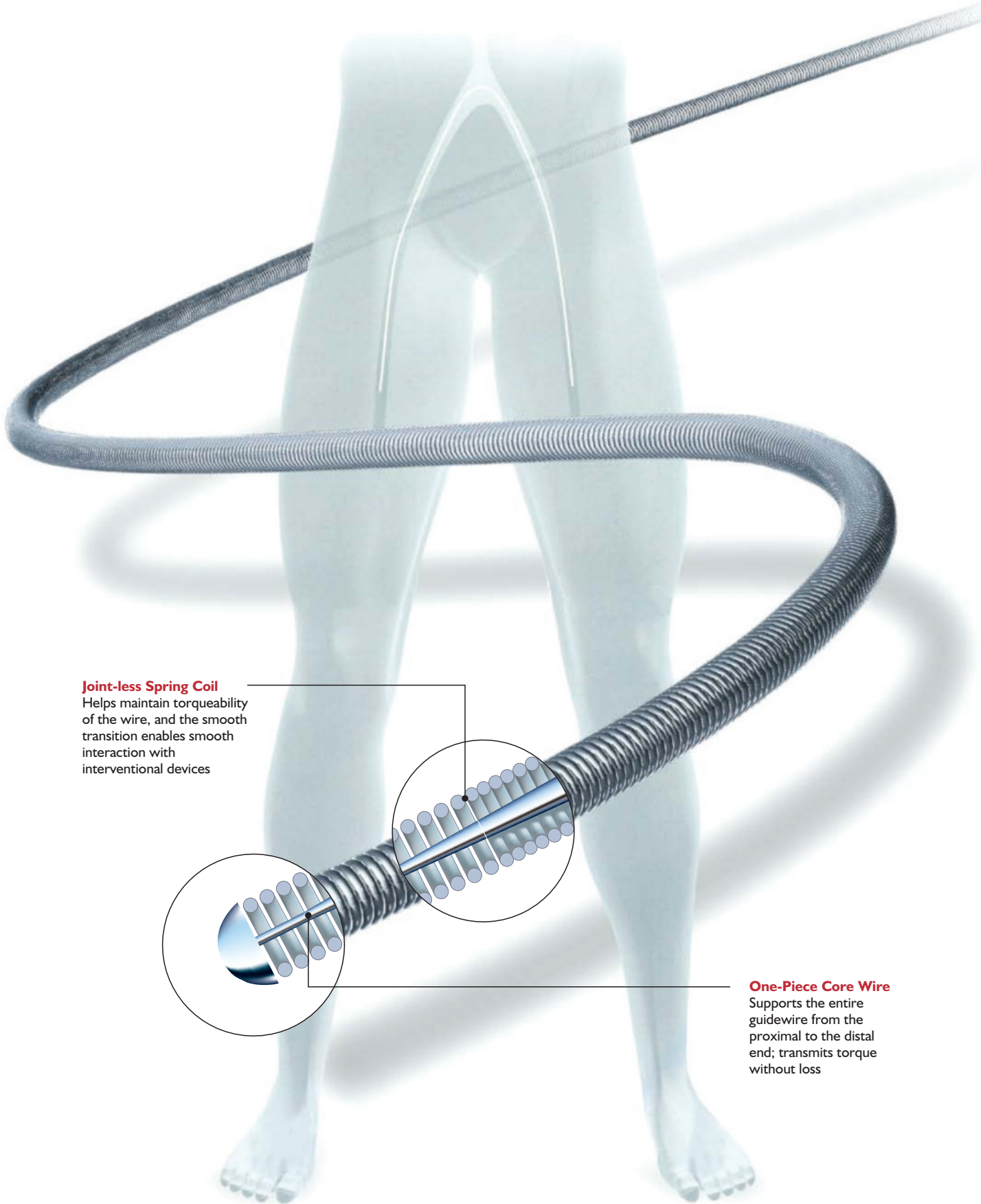

### Joint-less Spring Coil

Helps maintain torqueability of the wire, and the smooth transition enables smooth interaction with interventional devices

### One-Piece Core Wire

Supports the entire guidewire from the proximal to the distal end; transmits torque without loss

**Treasure Floppy 0.018"****Superior Handling and Distal Placement for Non-Occluded Lesions**

The Treasure Floppy has superior torque and handling characteristics with enhanced vessel trackability, which allows for improved distal placement and device delivery.

- Tip load – 4.0 g
- Tip radiopacity – 3 cm
- Total length – 190 cm, 300 cm
- Outer diameter – 0.018 inch
- SLIP COAT® Hydrophilic coating over the spring coil
- Spring coil length – 15cm
- PTFE coating over the shaft

**Treasure I2 0.018"****Controlled Torque and Manoeuvrability with Tactile Response**

The Treasure I2 is ideally suited as a CTO wire for the treatment of PAD, with exceptional control throughout the case due to the one-piece core and spring coil design.

- Tip load – 12.0 g
- Tip radiopacity – 15 cm
- Total length – 180 cm, 300 cm
- Outer diameter – 0.018 inch
- SLIP COAT® Hydrophilic coating over the spring coil
- Spring coil length – 15cm
- PTFE coating over the shaft

**Astato 30 0.018"****Superior Penetration Ability, Manoeuvrability with Tactile Response**

A high-penetration Hydrophilic guide wire with a tapered Hydrophobic tip and 30 gram tip load, tracks through calcified lesions whilst providing superb tactile feedback. This makes the Astato 30 the ideal choice for penetrating the proximal CTO cap.

- Tip load – 30.0 g
- Tip radiopacity – 15 cm
- Total length – 180 cm, 300 cm
- Outer diameter – 0.018 inch
- Tip outer diameter – 0.013 inch
- SLIP COAT® Hydrophilic coating over the spring coil
- Spring coil length – 15cm
- PTFE coating over the shaft

**Regalia XS 1.0 0.014"****Small Vessel Crossing with Advanced Slip Performance**

The Regalia XS 1.0 is a Polymer sleeve frontline wire for heavily stenosed, sub-total lesions and small BTK arteries. The SLIP COAT® Hydrophilic coating and Polymer sleeve provide slip performance and trackability, whilst maintaining excellent feedback and physician control. The Regalia XS 1.0 is also superb for tracking small and tortuous vessels.

- Tip load – 1.0 g
- Tip radiopacity – 3 cm
- Total length – 180 cm, 300 cm
- Outer diameter – 0.014 inch
- Polymer sleeve length – 22 cm
- Spring coil length – 12cm
- SLIP COAT® Hydrophilic coating over the polymer sleeve
- PTFE coating over the shaft


**Astato XS 20 0.014"****Superior Crossing Ability for Heavy Stenosis in Small Vessels**

Hydrophilic coated guide wire with a tapered Hydrophobic tip and 20 gram tip load. The Astato XS 20 increases the ability to cross complex lesions with heavy calcifications and / or tough, fibrous tissues, especially in smaller vessels.


- Tip load – 20.0 g
- Tip radiopacity – 17 cm
- Total length – 180 cm, 300 cm
- Outer diameter – 0.014 inch
- Tip outer diameter – 0.008 inch
- SLIP COAT® Hydrophilic coating over the spring coil
- Spring coil length – 17cm
- PTFE coating over the shaft

### Torque Performance

Asahi wires provide one-to-one torque performance for superior wire handling.

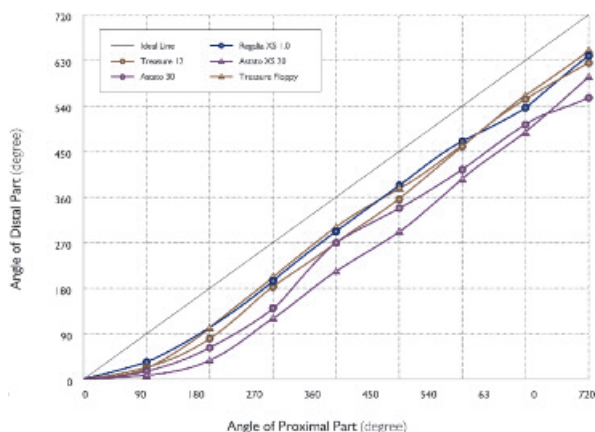

### Wire Selection Chart

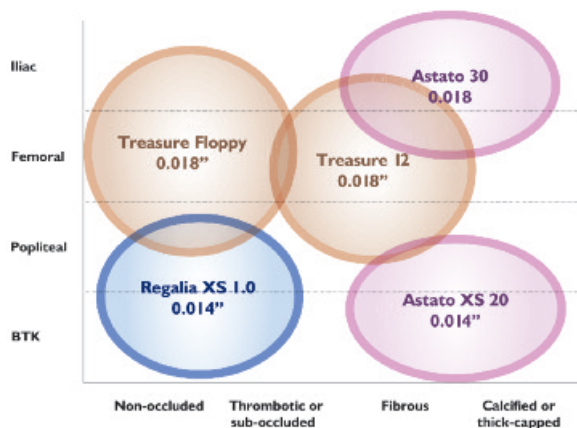

### Ordering information

Cat. Number	Diameter mm (inch)	Radiopaque Segment (cm)	Spring Coil Length (cm)	Length (cm)	Tip Shape
<b>Treasure Floppy</b>					
PAGHI8M072	0.45 (0.018)	3	15	190	Straight
PAGHI8M372	0.45 (0.018)	3	15	300	Straight
<b>Treasure 12</b>					
PAGHI8M070	0.45 (0.018)	15	15	180	Straight
PAGHI8M370	0.45 (0.018)	15	15	300	Straight
<b>Atrato 30</b>					
PAGHI8M071	0.45/0.33 (0.018/0.013)	15	15	180	Straight
PAGHI8M371	0.45/0.33 (0.018/0.013)	15	15	300	Straight
<b>Regalia XS 1.0</b>					
PAGPI40000	0.36 (0.014)	3	12	180	Straight
PAGPI40300	0.36 (0.014)	3	12	300	Straight
<b>Atrato XS 20</b>					
PAGHI43092	0.36/0.20 (0.014/0.008)	17	17	180	Straight
PAGHI43392	0.36/0.20 (0.014/0.008)	17	17	300	Straight

## Vascularperspectives

innovations for the interventionalist

Vascular Perspectives, Parkway House  
Palatine Road, Manchester, M22 4BD, UK

Tel: 0161 491 0585  
Fax: 0161 282 6540

Web: vascularperspectives.com  
Email: info@vpmed.co.uk